BBS

ES | EN | FR | PT | IT | DE

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ELECTRONIC BABY SCALE

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BALANÇA ELETRÓNICA DE BEBÉS

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V1.1
14/06/2022

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1. DEAR FRIENDS

Welcome to use our baby scale, the best growing companion! Please read the instruction carefully to ensure proper usage.

2. MAJOR SPECIFICATIONS & FEATURES


- Maximum capacity: 25kg/55lb; Division: 0.005kg/0.01 lb.
- Height Measurement range: 46-80cm/18.1-31.5in , division: 1 mm/0.1 in.
- Memory: 8 users supported, 15 records saved for each user. Data comparison mark in user mode.
- No.0 is guest mode. Data would not be record.

3. BATTERY INSTALLATION

Type	Diagram	Installation
4 X 1.5V AAA Batteries (Not included)		Install the AAA batteries properly by positive / negative pole

Note: Please use corresponding battery. Please remove battery if it is not in use in a long time (more than three months).

4. LCD SPECIFICATION

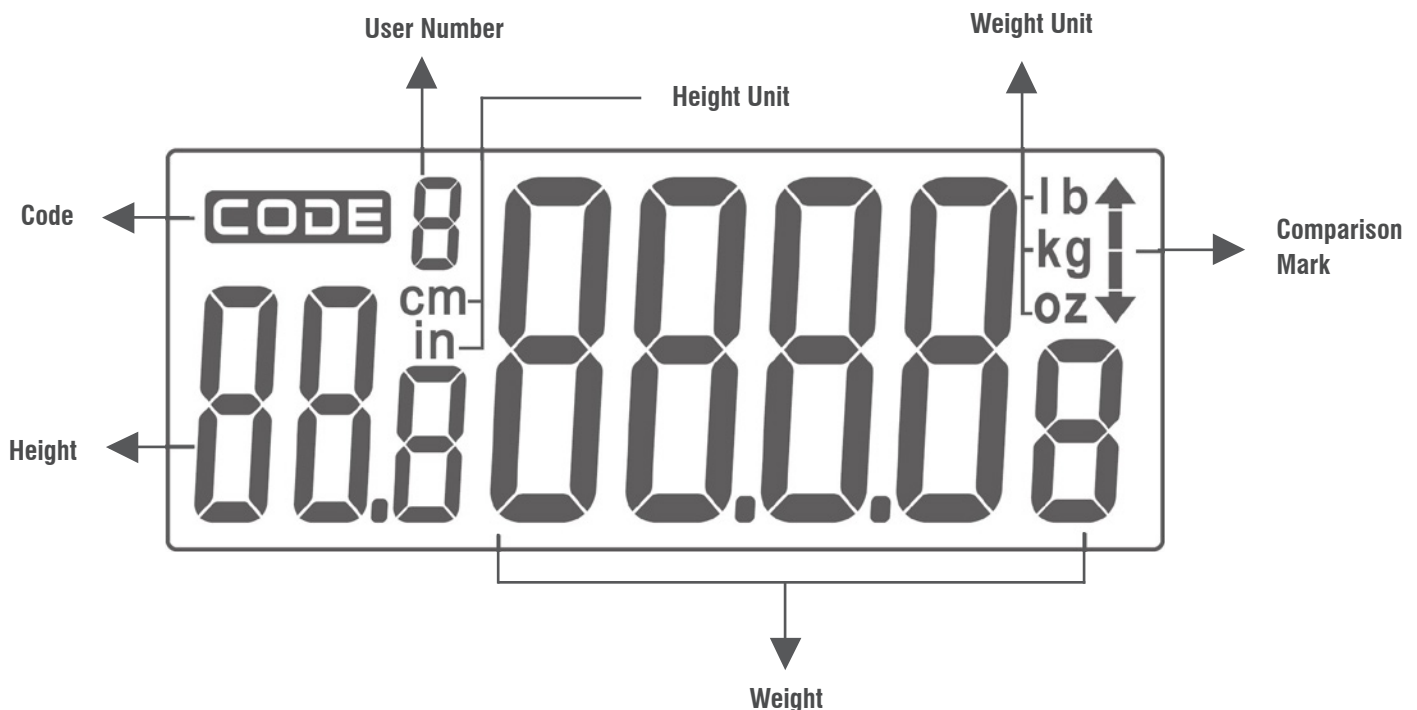5. BUTTON FUNCTIONS

Turn On/Off; ZERO.

Short press to hold / unhold measurements, long press to switch unit.

Short press access into history review, then go through every history record by every short press.

Long press access into user selection, then short press to choose the right user and confirm it by long pressing again.

6. HOW TO MEASURE WEIGHT AND HEIGHT

To ensure the accuracy of measurement, please use the scale on a flat and firm floor.

1. Press “” to switch on scale and go back to previous user number. Long press “” to choose user, and then short press to switch user, Long press again to confirm. Data will not be saved in guest mode, (user number “0”).

2. Pull out the silicone plug and insert the overhead board (Refer to Diagram 1 to 2).

3. Lift open up the handle to vertical for activating the height measurement, and keep it totally leaned toward initial position, LCD will show 46.0cm/18.1 in as initial value. If not, press “” to initialize the value, then pull out the measurement board to measure baby's height.

4. Lay down your baby carefully on the scale. The head should touch the overhead board, at the same time, the feet should touch the vertical handle.

5. Press “” to hold the measurements of weight and height. Results flash twice for confirmation and are stored in scale memory.



Comparison to previous data will be displayed in difference value repeatedly. In guest mode, measurements will be displayed for reference,

but will not be saved. (Press “” to cancel locking or let it turn off automatically without any operation.)

Note: About memory of scale: Maximum 8 user supported. 15 records for each user. The oldest record will be overwritten by the latest record when the records are more than 15.


Note: For weight measurement only, just skip step 2 & 3.

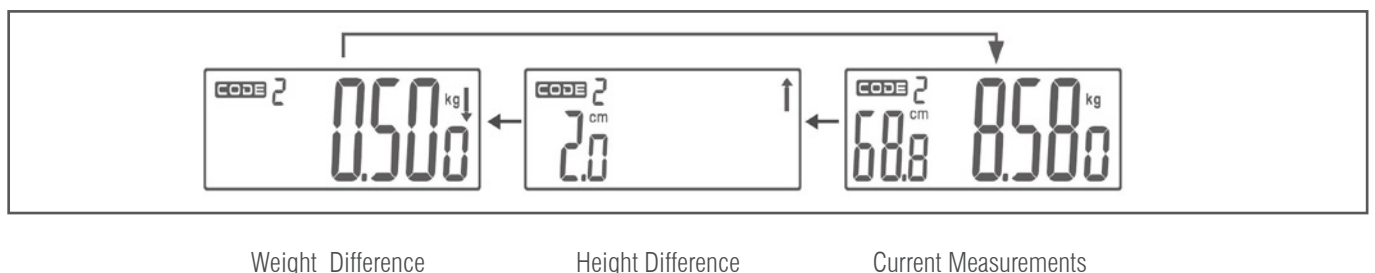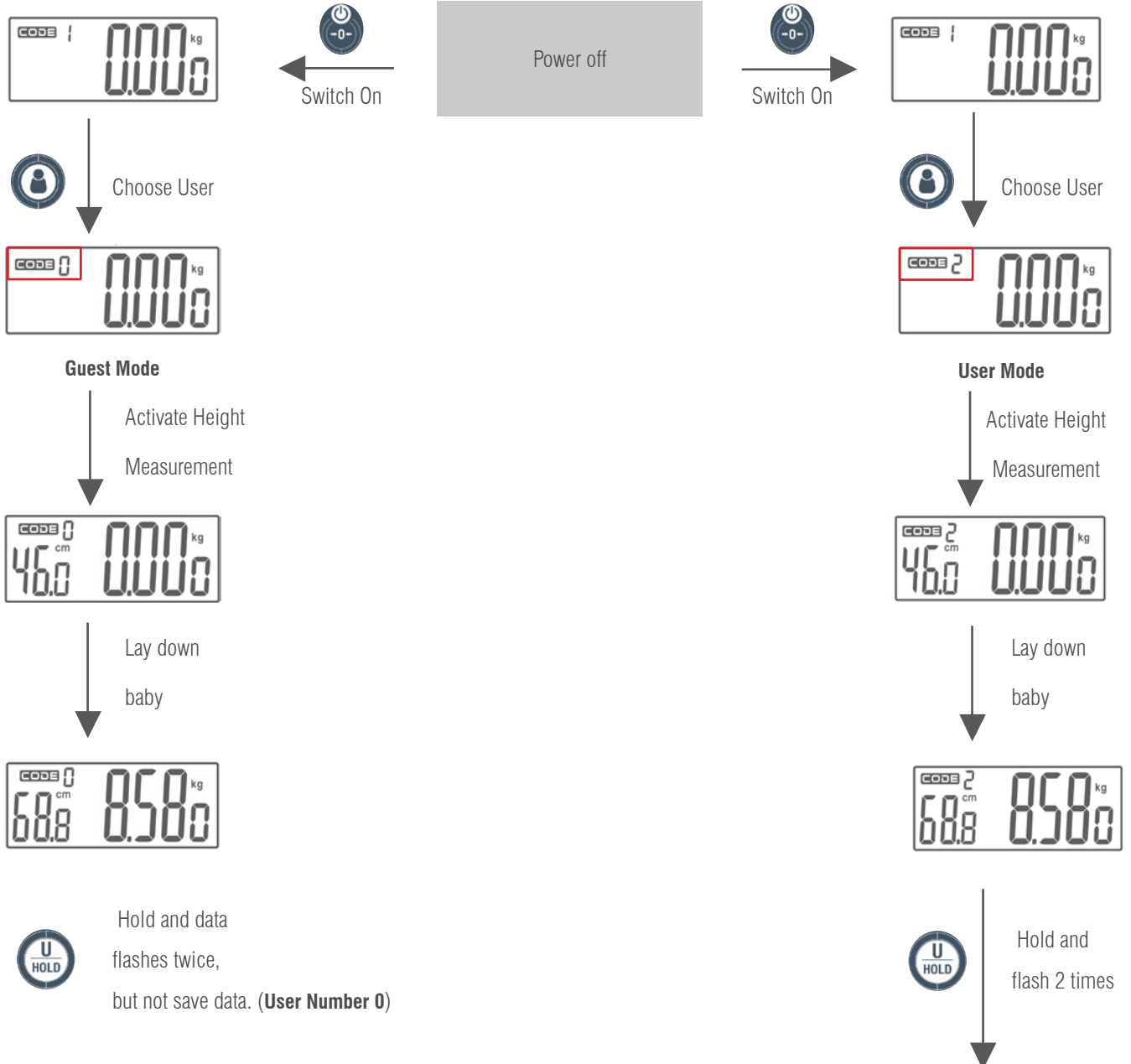Data saved, and repeatedly displayed (**User number different from 0**).

7. HISTORY

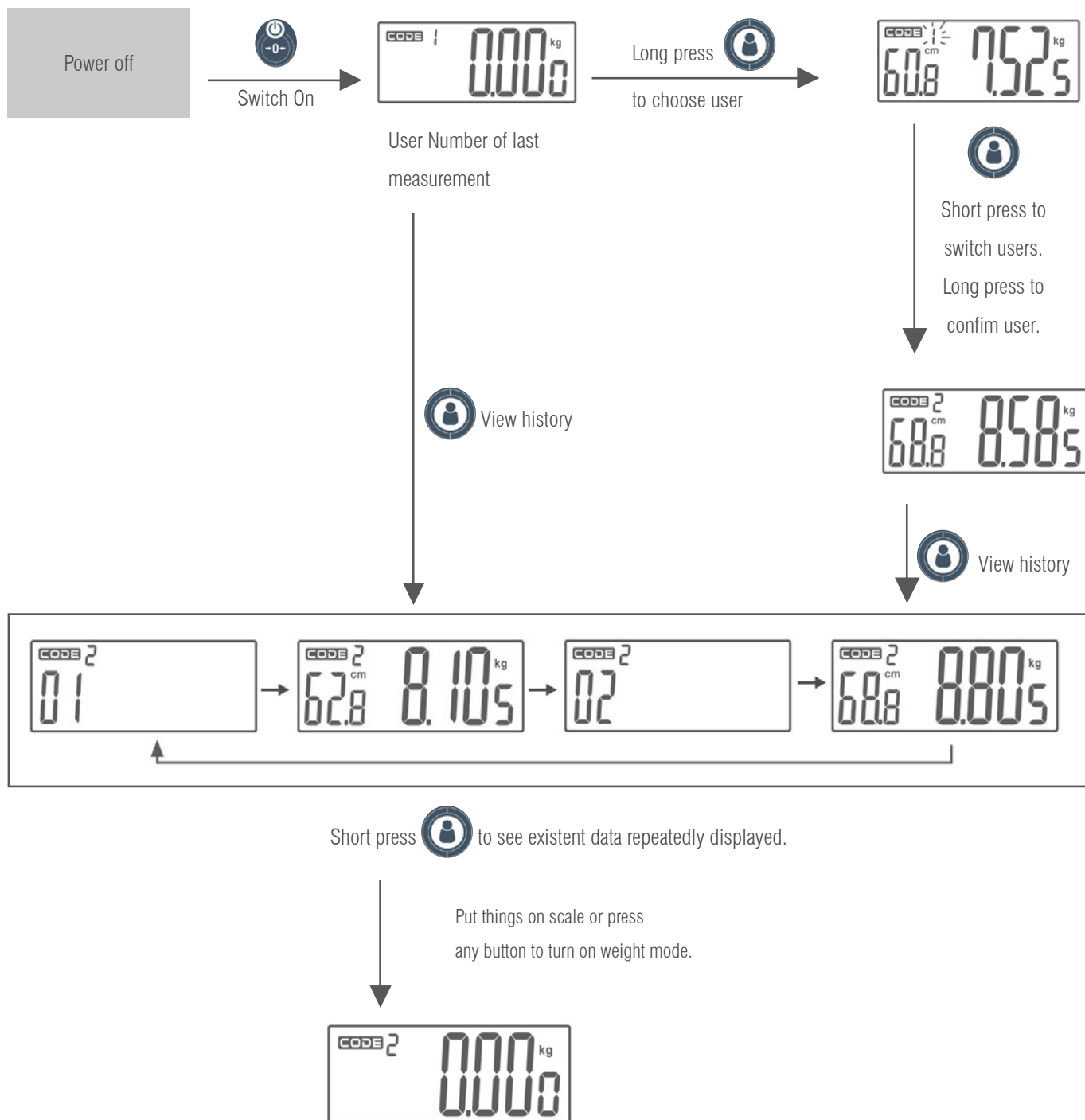Note: Long press “” to choose user.

Put things on scale or press any button to turn on weight mode.

8. AUTO OFF AND MANUAL OFF

For longer battery life, the scale will turn off automatically after about 1 minute without any operation, or long press “” to turn it off after use.

9. INDICATION

Warning indication	Name	Instructions and methods
Err	OVERLOAD INDICATION	Scale overload the max capacity (Max capacity is indicated on the scale.) Please slowly hold up to avoid damage.
Lo	LOW BATTERY INDICATION	Battery power is low. Please replace new battery.

10. CLEANING AND CARE

1. Clean the scale with a slightly damp cloth. DO NOT immerse the scale in water or use chemical/abrasive cleaning agents.
2. All plastics parts should be cleaned immediately if contact with fats, spices, vinegar and strongly flavored/colored foods, or acids citrus juices.
3. Always use the scale on a solid, flat surface.
4. Do not knock, impact and drop products.
5. This product is a precision instrument, please be careful.
6. Place in a cool dry place.
7. Check the battery if the scale can not be turned on.
8. Designed for family and non-commercial.

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trade mark propriety of | marca propiedad de | est une marque de | marca de propriedade de